


432Hertz Generator Micro-Current Technology

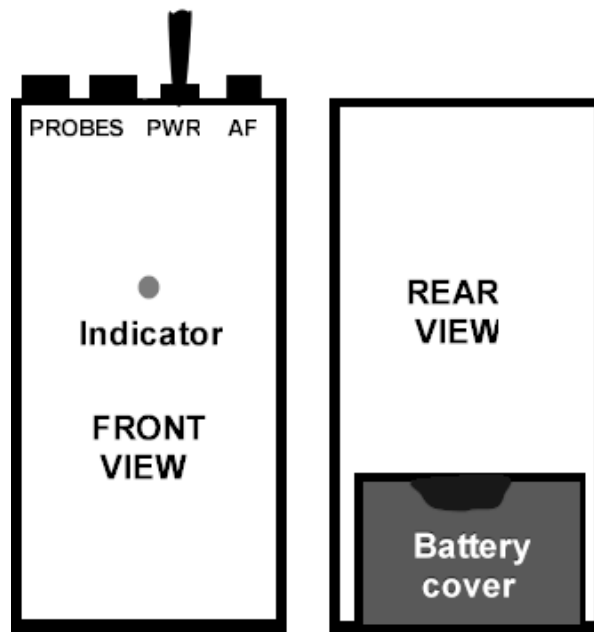

Specification:

432Hz modified micro-current waveform output +/- 0.03%

Size: 124 x 74 x 25

Weight: Generator - 150g including battery / Probes – 90g

Battery: 9 volt PM3/PP3 type

Probes: ~120 x 6mm copper

Probes Output: 432Hz :1.7VRMS <12mA@9v / 1Kohm DC

AF Output: @432Hz : 1VRMS/1Kohm

Operation:

Plug probes into RED and BLACK sockets on top of generator

Use switch located on top of generator labelled PWR to turn it ON or OFF.

The ON position will be indicated by a flashing indicator light on the front to the generator

The socket labelled AF is an output of 432Hz that can be connected to an amplifier for audible purposes by using a suitable cable. The AF socket requires a 3.5mm mono jack plug.

Change battery when indicator light appears dimmer than usual

To replace battery, slide battery cover (located on rear) downwards. Un-clip the battery and insert the new one. Observe the battery clips polarity.

Method #1

Hold probes in palms of hands. Turn generator ON

Maintain position for 30 minutes or as required

Method #2

With probes on the floor, place both feet over the probes (one probe under each foot).

Maintain position for 30 minutes or as required

Further useful information on 432 Hertz Healing Frequency, as well as Schumann Frequency can be found on the Internet.

NOTE: We manufacture equipment and are not practitioners. Use of the generator is entirely the responsibility of the user.